
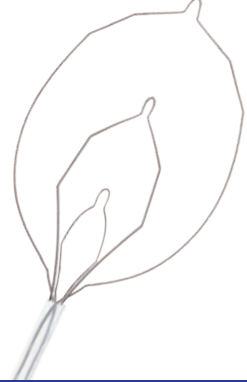

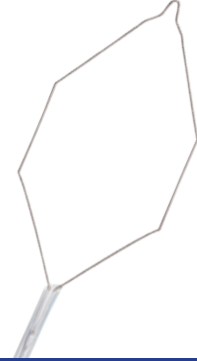
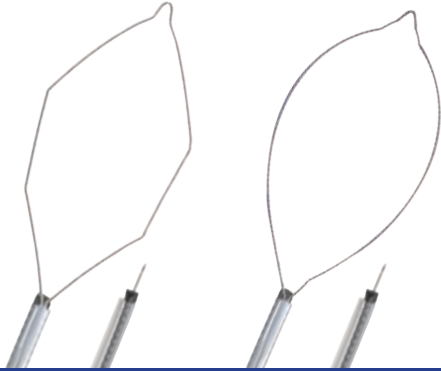

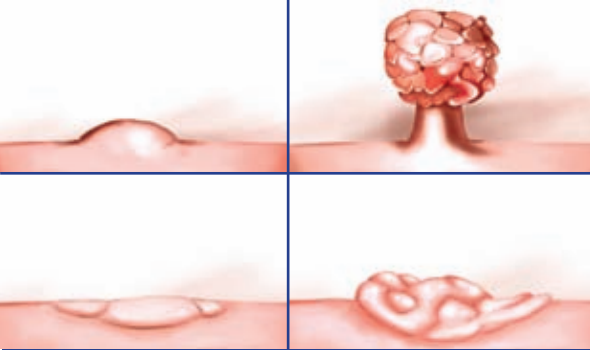





# US Endoscopy polypectomy portfolio reference guide

Snare/Device	Exacto® cold snare	Lariat® snare	Histolock™ resection device	Traxtion® snare	iSnare® system (oval & hex)
					
	BX00711115	BX00711119	BX00711117	BX00711113	BX00711085/ BX00711086 BX00711088/ BX00711089
Polyp Type	 <b>Diminutive (1-9mm)</b>	 <b>Multiple Polyps (all types)</b>	 <b>Very Flat/Sessile</b>	 <b>Large Flat/Sessile</b>	 <b>Large Flat/Sessile</b>
Snare Width	9mm	6mm, 10mm, 30mm	14mm	25mm	25mm
Snare Wire Diameter	Braided .012in. (.30mm)  33% thinner than standard snares (.018") for a clean cut	Braided .016in. (.41mm)	Monofilament flat wire	Braided .018in. (.46mm)	Braided .018in. (.46mm)
Features	<ul style="list-style-type: none"> <li>90% of polypectomies are less than 9mm in size<sup>1</sup>.</li> <li>Allows for full circumferential resection of polyp.</li> <li>Size and shape offer more control for placement around small polyps.</li> </ul>	<ul style="list-style-type: none"> <li>Three distinct shapes in one device.</li> <li>Offers versatility to handle cases with large polyps and multiple polyps of varying types and sizes.</li> </ul>	<ul style="list-style-type: none"> <li>Monofilament stiff wire easily recruits and resects difficult flat polyps.</li> <li>Ideal for challenging EMR cases with scarred or adherent polypoid tissue.</li> </ul>	<ul style="list-style-type: none"> <li>Braided stiff wire that grips mucosa and maintains grip on larger flat polyps.</li> <li>Ideal for en bloc resection of large flat lesions.</li> </ul>	<ul style="list-style-type: none"> <li>Injection needle and snare in one device operated by a multi-functional handle.</li> <li>Ideal for saline assisted polypectomy cases.</li> </ul>

<sup>1</sup>Tolliver KA, Rex DK. Colonoscopic Polypectomy. Gastroenterol Clin N Am 37. 2008; 229-251.

Refer to the current **Instructions For Use** for complete product usage instructions.

US Endoscopy is a subsidiary of STERIS Corporation.  
© US Endoscopy or its affiliates. All rights reserved.  
Listening...and delivering solutions®, the US Endoscopy listening...and delivering solutions design, and all marks denoted with ® or ™ are registered with the U.S. Patent and Trademark Office, or trademarks owned by US Endoscopy.